

Cooled incubators (ST)

Application

- BOD determination
- microbiological research
- plant growing and microorganisms breeding at specified temperature
- storage of liquids and samples for physicochemical analysis



Excelsior Scientific



Cooled incubators (ST) can provide a stable temperature between +3...+70°C regardless of ambient conditions.

Calibration



All thermostatic equipment manufactured by Excelsior Scientific can be provided with Calibration Certificate issued by our accredited Measurement Laboratory. Detailed information on accreditation can be found on our website: www.excelsiorscientific.com






The **BASIC, COMFORT, PREMIUM** models are equipped with a PID microprocessor controller with an LCD graphic display and illuminated touch buttons.

Controller advantages


- six segment temperature-time profile
- loop function up to 99 times or endless
- 3 user programs memory
- adjustable start delay feature (from 1 min to 99:59 h)
- adjustable hold at set point time for temperature and lighting (for ST/FOT) from 1 min to 31 days / 1 min to 99:59 h or continuous operating
- recording of min, average and max temperature value for each segment
- overview of set and current parameters while operating
- audible and visual temperature alarm
- operating with temperature priority mode
- defrosting function
- temperature sensor fail alarm
- power failure control system (program continued after restoring power)
- digital timer
- real time clock
- auto-diagnostic function
- internal memory to store up to 2046 data records
- forced air convection with optional fan speed control 50-100%
- automatic fan shut-down after completing the program

Detailed description of parameters on page 82.

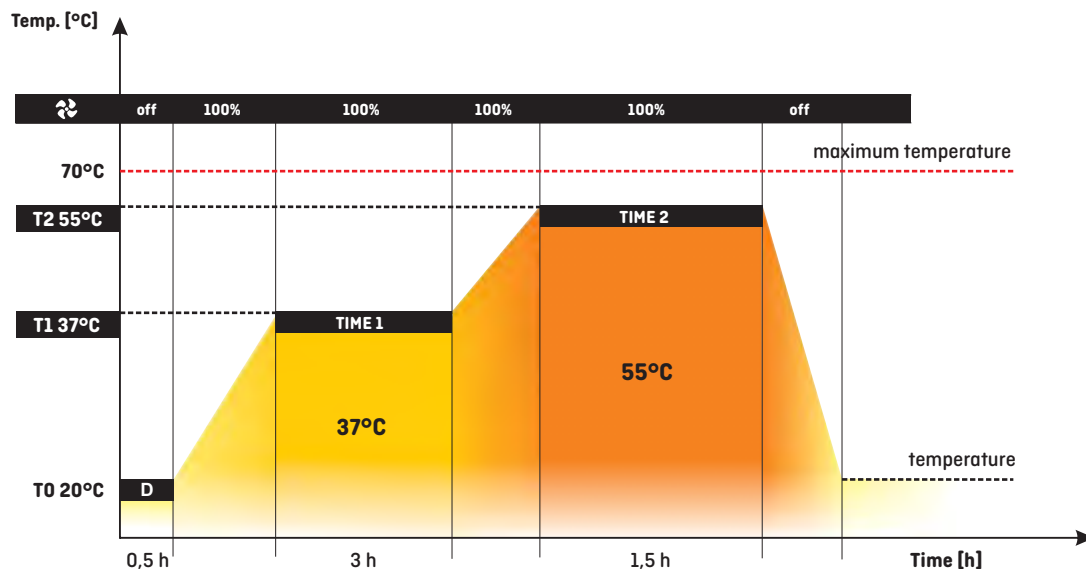

parameter set by the User D - adjustable start delay (1 min...99:59 h) - fan speed control (option) T1/T2/ - segment temperature TIME 1/2/3 - segment time

Control panel

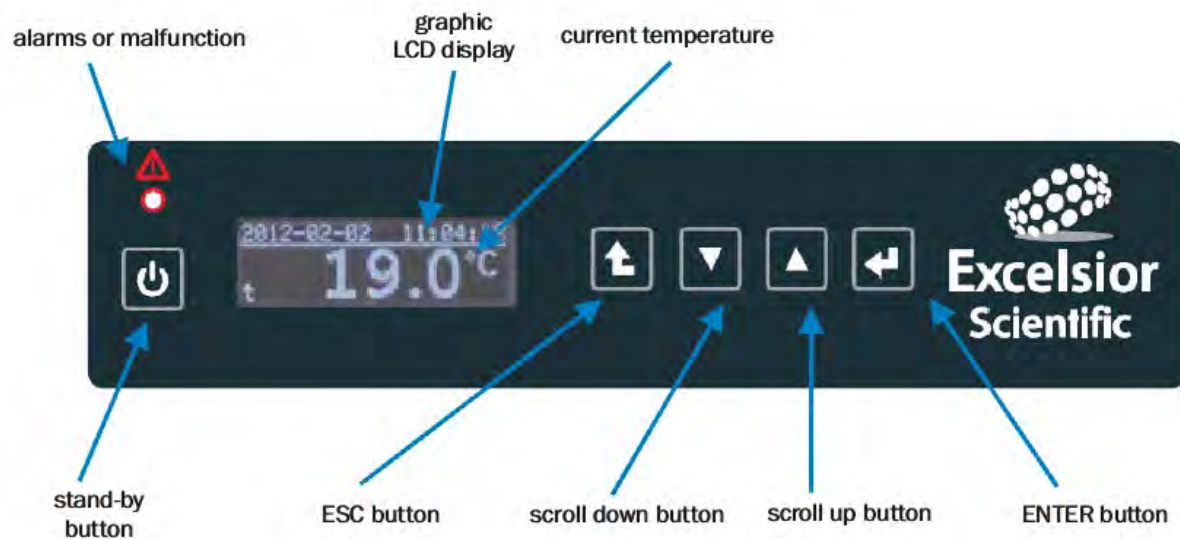

Standard features

- temperature range +3...+40°C
- quality control protocol (at +37°C)
- English instruction manual
- available menu languages: Czech, English, Estonian, French, German, Italian, Latvian, Polish, Portuguese, Russian, Spanish
- temperature protection 1.0 class for BASIC and COMFORT models and 2.0 for PREMIUM models to DIN 12880
- open door alarm
- wheels in standard for models ST 1200 and 1450

RS 232 and USB ports for data transfer

internal LED light

access port: Ø30 mm

door lock

wire shelves for BASIC and stainless steel wire shelves for COMFORT and PREMIUM models

solid door



PREMIUM TOP+ version


All the units in TOP+ version are equipped with PID microprocessor controller with a large (5,7") full colour touch screen, intuitive menu and user friendly software. They can be connected to Ethernet network for remote control from any computer, being one of the greatest advantages.

Controller advantages


- multi-segment temperature-time profile (up to 100)
- loop function up to 99 times or endless
- adjustable start delay feature (from 1 min to 99:59 h or date/time)
- access control via login
- 7-days programming
- adjustable hold at set point time for temperature and lighting (for ST/FIT) from 1 min to 999:59 h, or continuous operating
- adjustable ramps
- overview of set and current parameters while operating
- recording of min, average and max temperature value for each segment
- Administrator function to manage User accounts
- possibility of temperature calibration by the User
- audible and visual temperature alarm
- operating in temperature or time priority mode
- defrosting function
- temperature sensor fail alarm
- power failure control system (program continued after restoring power)
- digital timer
- real-time clock
- auto-diagnostic function
- forced air convection with fan speed control (50-100%)
- automatic fan shut-down after completing the program

GLP supporting functions

- password protected settings
- 20 user programs memory
- internal memory to store up to 4100 data records for each User, possibility to overview the values on the display or a PC computer in tabular or graphic form
- USB port to allow direct data recording or transfer into a flash drive
- events registry

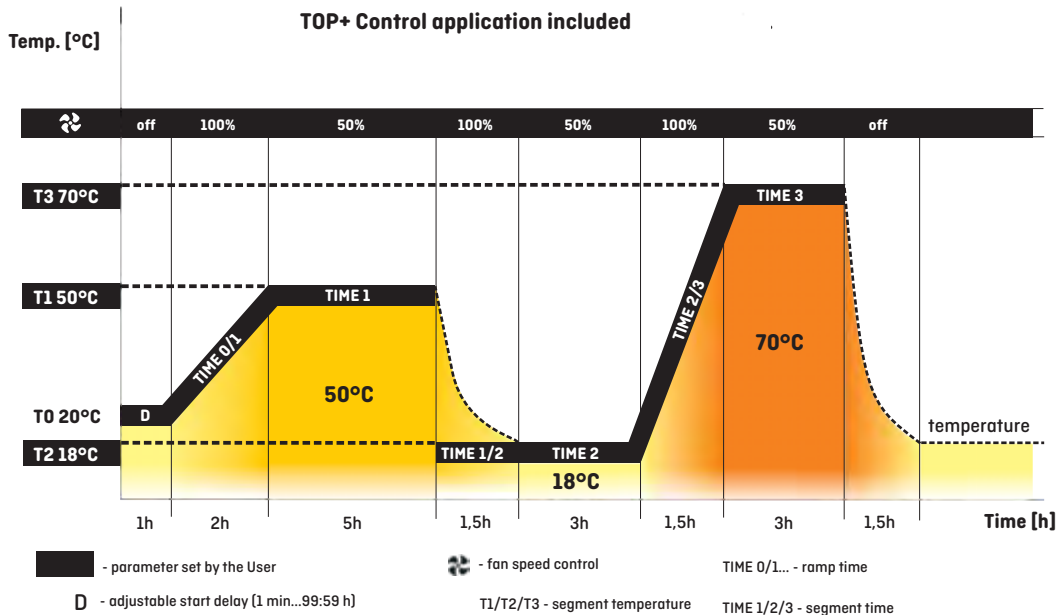

Control panel

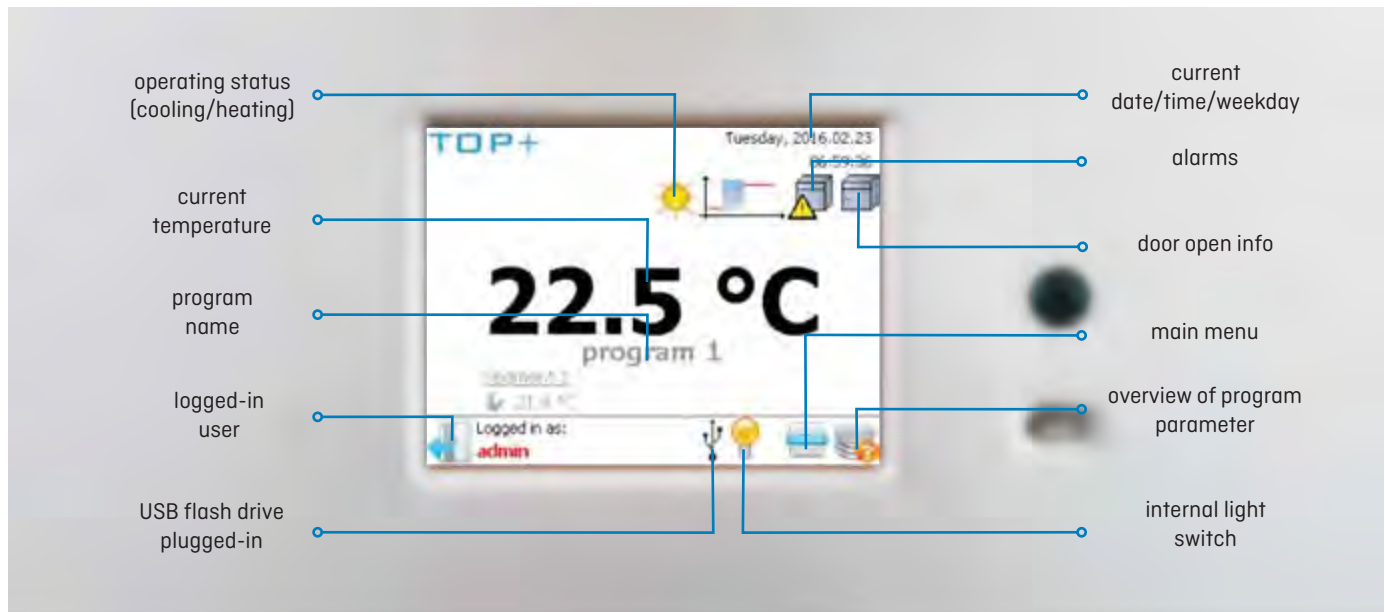

Standard features

- temperature range +3...+70°C
- Ethernet cable
- TOP+ Control software
- quality control protocol (at +37°C)
- English instruction manual
- available menu languages: English, Estonian, French, German, Hungarian, Italian, Latvian, Polish, Portuguese, Romanian, Russian, Spanish
- temperature protection 3.3 class to DIN 12880
- open door alarm
- wheels in standard for models ST 1200 and 1450

door lock

USB port to allow direct recording and data transfer onto a flash drive

RS 232 interface and LAN port

internal LED light

access port: Ø30 mm

wire stainless steel shelves

solid door



Cooled incubators (ST)

ST

		ST 1	ST 2	ST 3	ST 4	ST 5	ST 6	ST 500	ST 700	ST 1200	ST 1450
Parameter											
air convection		forced									
chamber capacity [l]		70	150	200	250	300	400	500	625	1365	1460
working capacity [l]		55	122	163	203	243	324	386	450	1229	1307
door type		solid / glass or double ¹ (option)									
temperature range [°C]		+3...+40 / up to +70 (option) / +3...+70 in PREM TOP+									
temperature resolution [°C]		every 0,1									
controller		microprocessor with external LCD graphic display									
interior	BASIC	aluminum									
	COMF	stainless steel to DIN 1.4016									
	COMF/S	stainless steel to DIN 1.4016									
	PREM (TOP+)	stainless steel to DIN 1.4301									
	PREM/S (TOP+)	stainless steel to DIN 1.4301									
housing	BASIC	powder coated sheet									
	COMF	powder coated sheet									
	COMF/S	polished stainless steel									
	PREM (TOP+)	powder coated sheet									
	PREM/S (TOP+)	polished stainless steel									
overall dims ² [mm]	A width	570	620	620	620	620	620	660	750	1480	1450
	B height	600	860	1060	1260	1460	1860	1990	1990	1990	1970
	C depth	680	650	650	650	650	650	810	860	860	950
internal dims ³ [mm]	D width	430	480	480	480	480	480	430	480	2x480	2x490
	D' width	470	520	520	520	520	520	510	600	1310	1340
	E height	430	660	860	1060	1260	1660	1510	1510	1510	1460
	F depth	300	420	420	420	420	420	650	690	690	750
	F' depth	360	480	480	480	480	480	-	-	-	-
	G depth	-	320	320	320	320	320	-	-	-	-
	H height	-	440	640	840	1050	1440	-	-	-	-
	I height	-	-	-	-	-	-	1380	1360	1360	1300
max shelf workload ⁴ [kg]	-	10	10	10	10	10	10	20	30	30	30
	Pw ⁵ version	on request						100	100	100	100
max unit workload [kg]	-	20	30	40	50	60	60	100	150	300	300
	W ⁶ version	on request									
nominal power [W]		160	170	170	330	330	330	400	400	550	550
weight ⁷ [kg]		32	54	59	69	75	90	105	115	185	200
over temperature protection		class 1.0 to DIN 12880 / class 3.3 (option) / class 3.3 in PREM TOP+									
power supply*		230V 50 Hz									
shelves fitted/max		2/2	3/4	3/4	4/6	4/7	4/10	3/11	3/11	2 x 3/11 ⁸	2 x 3/11 ⁸
warranty		24 months									
manufacturer		Bioscience Resource Project									

all the above technical data refer to standard units (without optional accessories)

* - 230V 60Hz, 115V 60Hz also available

1 - additional internal glass door

2 - ST 1-6 in TOP+ version are 60 mm higher, depth doesn't include 50 mm of power cable

3 - dims of units with double door can be smaller

4 - on uniformly loaded surface

5 - reinforced shelf







6 - reinforced version

7 - for units in BASIC version with solid door

8 - two columns with 3 shelves each

Options and accessories



		ST 1/1	ST 1/1/1	ST 2/2	ST 2/3	ST 2/4	ST 3/3	
Parameter								
air convection		forced						
chamber capacity [l]		70/70	70/70/70	150/150	150/200	150/250	200/200	
working capacity [l]		55/55	55/55/55	122/122	122/163	122/203	163/163	
door type		solid / glass or double ¹ (option)						
temperature range [°C]		+3...+40 / up to +70 (option) / +3...+70 in PREM TOP+						
temperature resolution [°C]		every 0,1\$						
controller		microprocessor with external LCD graphic display						
interior		aluminum						
		BASIC						
		COMF						
		COMF/S						
		PREM (TOP+)						
housing		stainless steel to DIN 1.4016						
		stainless steel to DIN 1.4016						
		stainless steel to DIN 1.4301						
		stainless steel to DIN 1.4301						
		stainless steel to DIN 1.4301						
overall dims ² [mm]		powder coated sheet						
		powder coated sheet						
		polished stainless steel						
		powder coated sheet						
		polished stainless steel						
internal dims ³ [mm]		A width	570	570	620	620	620	620
		B height	1170	1740	1680	1880	2080	2080
		C depth	680	680	650	650	650	650
		D width	470	470	520	520	520	520
		D' width	470	470	520	520/520	520/520	520
		E height	430	430	660	660/860	660/1060	860
		F depth	300	300	420	420	420	420
		F' depth	360	480	480	480/480	480/480	480
		G depth	-	320	320	320	320	320
max shelf workload ⁴ [kg]		H height	-	440	640	820	1050	1440
		-	10	10	10	10	10	10
		Pw ⁵ version	on request					
max unit workload [kg]		-	20	20	30	30/40	30/50	40
		W ⁶ version	on request					
nominal power [W]		320	480	350	350	350	350	
weight ⁷ [kg]		65	98	109	114	124	119	
over temperature protection		class 1.0 to DIN 12880 / class 3.3 (option) / class 3.3 in PREM TOP+						
power supply*		230V 50 Hz						
shelves fitted/max		see table for single chamber models						
warranty		24 months						
manufacturer		8&VXf beFVXagVV						

all the above technical data refer to standard units (without optional accessories)

* - 230V 60Hz, 115V 60Hz also available

1 - additional internal glass door

2 - depth doesn't include 50 mm of power cable

3 - dims of units with double door can be smaller


4 - on uniformly loaded surface

5 - reinforced shelf

6 - reinforced version

7 - for units in BASIC version with solid door


Options and accessories

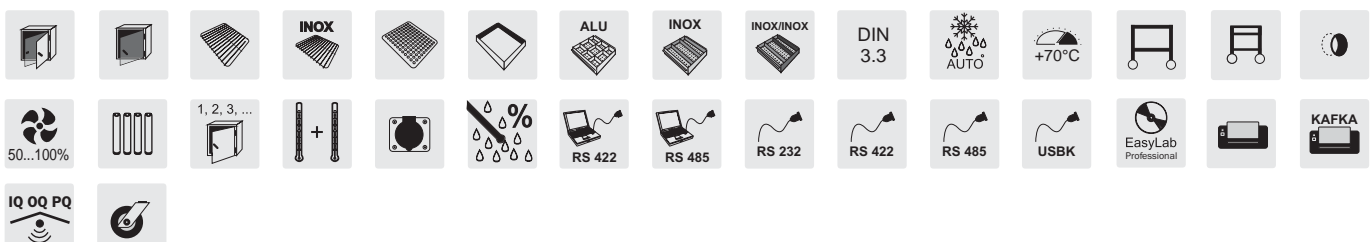

Cooled incubators (ST)


ST

ST 1

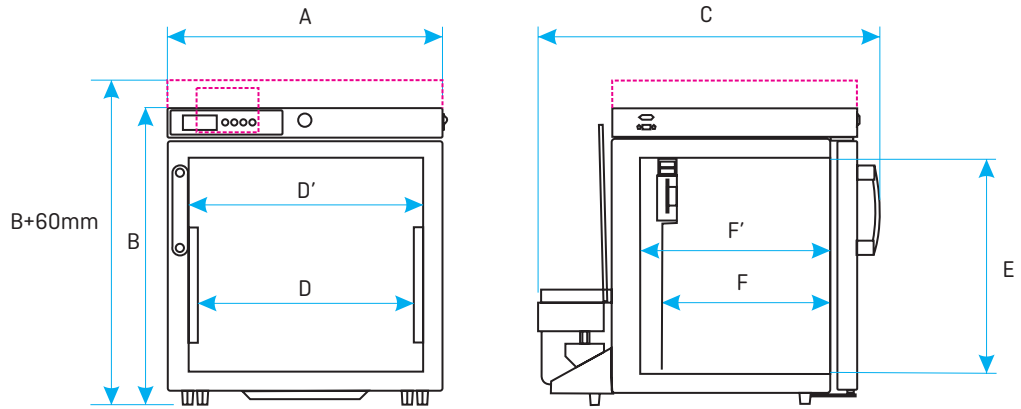


----- - TOP+ version

ST 2/3/4/5/6

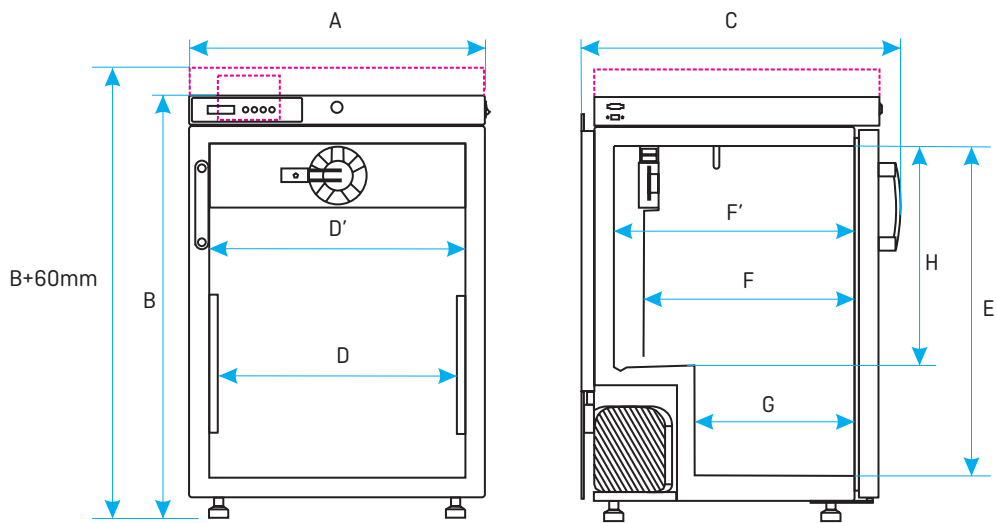

----- - TOP+ version

ST 500/700/1200/1450

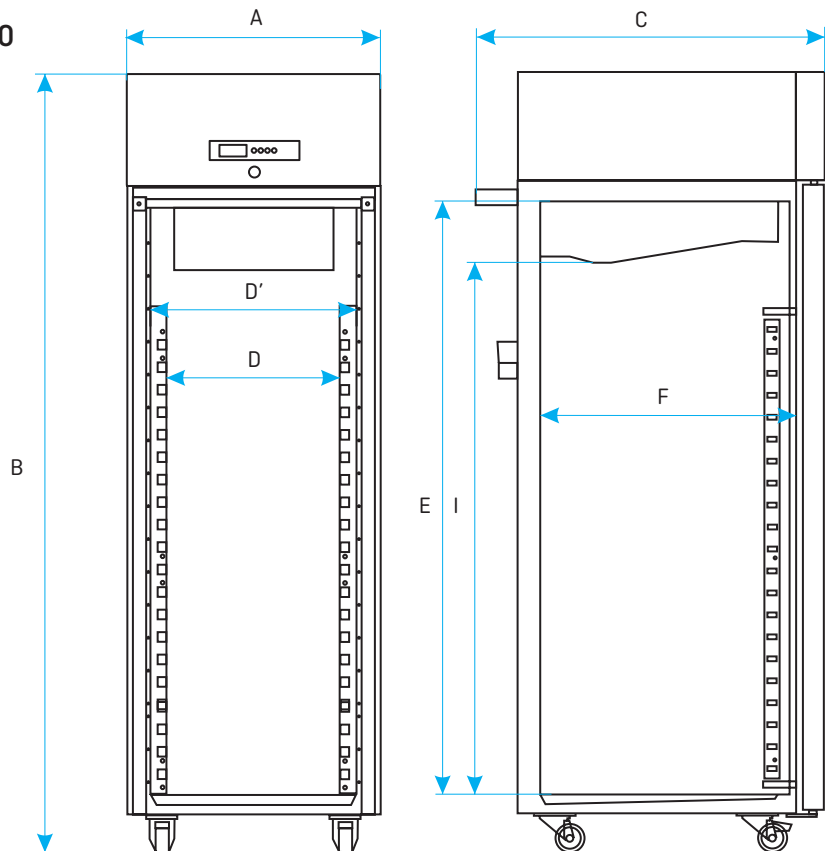

Cooled incubators (ST) with photoperiodic system

▲ Cooled incubators (ST) with photoperiodic system

The photoperiodic (FOT) and phytotron (FIT) systems allow day and night simulation. The FOT option allows the lights to be turned on or off, while the FIT option features additional intensity control.

The BASIC, COMFORT and PREMIUM versions of cooled incubators (ST) can be equipped with the FOT option, while the PREMIUM TOP+ version (ST 500, 700, 1200, 1450) with the FIT system.


▲ Program possibilities with FOT option

- day and night simulation software to control light (on/off), time and temperature separately for each segment
- temperature range for "night" simulation: +3 ... +50°C
- temperature range for "day" simulation: +10 ... +50°C
- lamps installed on side walls
- fluorescent lamp 840 type (daylight) as standard
- operating with time priority (see page 82)

▲ Photoperiodic system (* /FOT option) for single and double chamber cooled incubators (ST)**

Option	ST/FOT2	ST/FOT4	ST/FOT6	ST/FOT8	ST/FOT10	ST/FOT15
available for models	ST1 ST 1/1	ST 2 ST 2/2	ST 2; ST 3 ST 2/2 ST 3/3	ST 4 ST 5	ST 500 ST 700	ST 1200 ST 1450
temperature range with photoperiod ON [°C]	+10 ... +50					
number of lamps on walls	2	4	6	8	10	15 (3 columns with 5 pieces)
adjustable illumination intensity	no	no	no	no	no	no

** for ST models with */FOT option, inner dims can be narrower by 4 cm on each side. FOT option is factory preinstalled. There is no possibility to order it separately.

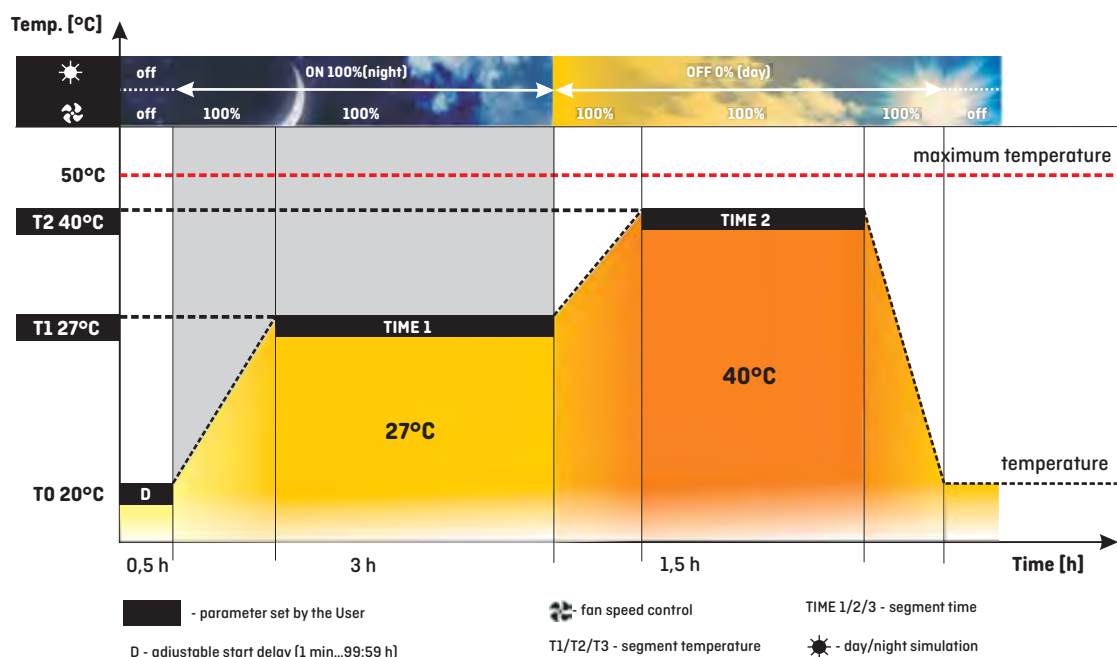

Cooled incubators (ST) with phytotron system

The PREMIUM TOP+ version of cooled incubators (ST 500, 700, 1200, 1450) can be equipped with the FIT system.


Program possibilities with FIT option

- day and night simulation software to control light intensity [%], time temperature and fan speed separately for each segment
- temperature range for „night” simulation: +3...+60°C
- temperature range for „day” simulation: +10... +50°C
- lamps installed in over-shelf panel (FIT P), side walls (FIT S), door (FIT D), door and side walls (FIT DS)
- fluorescent lamp 840 type (daylight) used as standard
- temperature or time priority program (see page 82)

Phytotron system (* /FIT option) for cooled incubators (ST) (PREMIUM TOP+)

Option**	ST/500/700/FIT DS	ST/500/700/FIT P	ST/500/700/FIT S	ST/1200/FIT P	ST/1450/FIT P
temperature range with photoperiod ON [°C]	+10 ... +50°C				
number of over-shelf panels with illumination std/max	-	1/3	-	1/3	1/3
lamps in door	yes	no	no	no	no
lamps in walls	yes	no	yes	no	no
adjustable illumination intensity	yes	yes	yes	yes	yes

** FIT DS - lamps in door and walls; FIT S - lamps in walls; FIT P - over-shelf panels

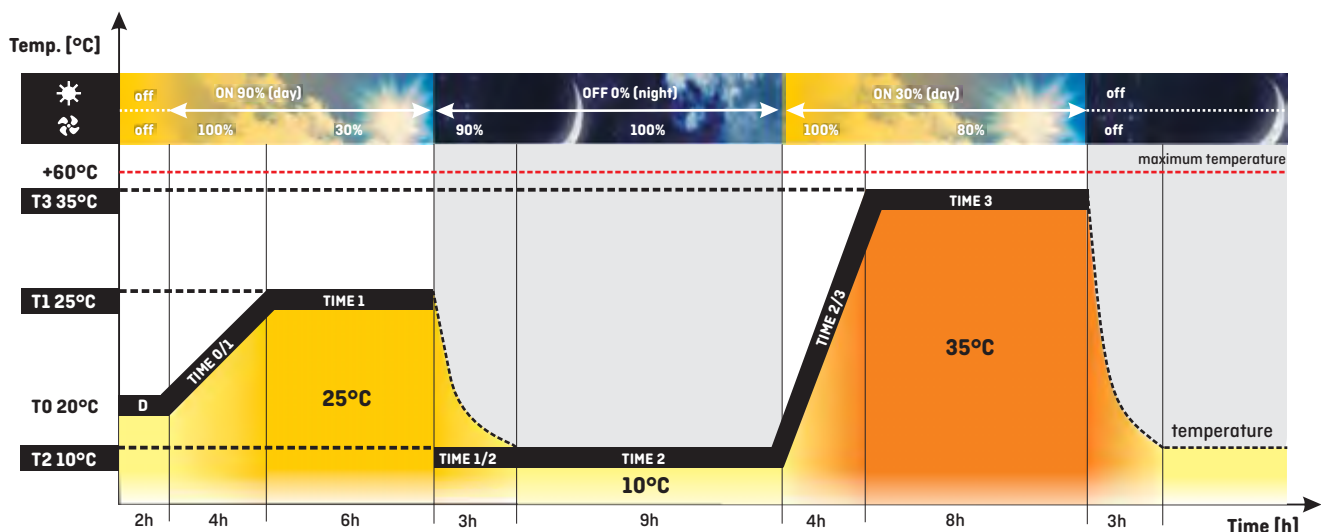

- - parameter set by the User
- ☀ - day/night simulation with illumination intensity control
- ⌚ - adjustable start delay (1 min...99:59 h)
- ⏸ - fan speed control
- ⏱ - TIME 0/1... - ramp time
- ⏱ - TIME 1/2/3 - segment time
- ⏱ - T1/T2/T3 - segment temperature